

CELLER[®]
Max

USER MANUAL



ACCELERATE YOUR WELLNESS POTENTIAL

WWW.CELLER8.COM

Contents

Introduction	3
Warnings & Safety Information	4-5
Contraindications	6-7
Precautions	8
Package Contents	9
Using the CELLER8 [®] Max	10-13
Remote Control Function	14-15
Applying the Applicators	16 -17
Technical Specifications	17
Maintenance & Cleaning	18-19
Troubleshooting	19-20
Warranty & Support	21
Customer Support Details	22

Introduction

Thank you for choosing the CELLER8® Max, a next-generation high-intensity Pulsed Electromagnetic Field (PEMF) device. Designed for personal and professional use, the CELLER8® Max delivers powerful PEMF sessions intended to support overall wellbeing, energy, and recovery through the application of targeted electromagnetic pulses.

PEMF technology has been explored for decades and is widely used in wellness, fitness, and performance recovery environments around the world. The CELLER8® Max builds on this tradition by offering a robust, user-friendly platform for those seeking a high-powered PEMF experience.

Important Notice:

The CELLER8® Max is not a medical device. It is not intended to diagnose, treat, cure, or prevent any medical condition or disease. No medical claims are made or implied. Always consult with a qualified healthcare provider before starting any new wellness or recovery protocol, especially if you have a pre-existing condition or health concern.

This manual will guide you through important safety information, operating instructions, and best practices to help you get the most out of your CELLER8® Max sessions.

Warnings & Safety Information

Before using the CELLER8[®] Max, please carefully review the following safety guidelines. These precautions are intended to help ensure safe and effective operation of the device and to prevent any risk of injury or equipment damage.

General Safety Guidelines

Read through the manual thoroughly before first use. Familiarise yourself with all controls, connections, and attachments.

Use only as intended. The CELLER8[®] Max is designed exclusively for personal and professional wellness use. It is not a substitute for medical advice or treatment.

Do not attempt to modify or repair the device. Any unauthorised servicing may void your warranty and compromise safety.

Do not operate the device if it appears damaged or if the power cable is frayed or compromised.

Supervise children or vulnerable individuals during use. The device is not intended for unsupervised use by anyone under 18 years of age.

Please leave at least 20 minutes for the device to cool, between sessions or users. This will avoid overheating of the applicators and the internal components of the device. Misuse of this cooling period will void any remaining warranty on both the device, fans or applicators.

Unplug the device when not in use, and store in a dry, cool, dust-free environment.

Please clean the device between users, using the cleaning guidelines later in this manual.

Do not use without supervision of a trained operator

Electromagnetic & Electrical Safety

Do not use near medical diagnostic or monitoring equipment. This includes ECG monitors, EEG machines, or similar devices.

Keep a minimum distance of 2 meters between the CELLER8[®] Max and other electronic devices, including mobile phones, laptops, TVs, and Wi-Fi routers. High-intensity electromagnetic fields can interfere with sensitive electronics.

Do not use with extension cords or power strips unless they are surge-protected and rated for high-powered electronics.

Remove all cards, wallets, purses, and electronic items from pockets and keep them away from the device before use, as the device may interfere with or cause damage to these items.

Environmental Conditions

Indoor use only. The device must be operated in a clean, dry, well-ventilated environment.

Do not expose to water, high humidity, or direct sunlight.

Avoid placing the device on soft surfaces (e.g., beds, couches) that may block ventilation or airflow.

During Use

Stop use immediately if you experience dizziness, nausea, palpitations, or any discomfort.

Do not sleep while the device is operating.

Do not use the device while driving or operating machinery.

Avoid contact with metal objects such as zippers, jewellery, or watches during sessions, as these may conduct energy and cause localised heating or discomfort.

Contraindications

The CELLER8® Max delivers high-intensity pulsed electromagnetic fields, which are not suitable for all individuals. Use of the device is contraindicated in the following cases:

Implanted Electronic Medical Devices

Do not use the CELLER8® Max if you have any type of implanted electronic device, including but not limited to:

- Pacemakers
- Cochlear implants
- Insulin pumps
- Neurostimulators or spinal cord stimulators
- PEMF can interfere with the function of these devices, potentially leading to serious health risks or device malfunction.

Pregnancy

The CELLER8® Max should not be used by individuals who are pregnant or suspect they may be pregnant. The effects of high-intensity PEMF on fetal development have not been studied sufficiently, and safety cannot be guaranteed.

Epilepsy or Seizure Disorders

PEMF exposure may affect neurological activity. Individuals with a history of epilepsy, seizures, or unexplained loss of consciousness should avoid using this device, as a precaution.

Active Bleeding or Hemorrhagic Conditions

If you have a condition that involves active internal or external bleeding, such as gastrointestinal bleeding, open wounds, or recent surgery, PEMF use may increase circulation and exacerbate bleeding. Consult your healthcare provider before use.

Organ Transplants

If you have received an organ transplant, PEMF may influence immune modulation and tissue behaviour in ways that are not compatible with anti-rejection medication. Always consult your transplant team before considering PEMF therapy.

Severe Cardiac Conditions

Individuals with serious heart conditions, such as arrhythmias, heart failure, or those who have had recent cardiac surgery, should avoid high-intensity PEMF use unless cleared by a cardiologist. While PEMF does not directly target the heart, electromagnetic fields may affect electrical signals in the body.

Cancer or Tumours

If you have a known or suspected malignancy (cancerous growths), do not use the device unless you are under the supervision and recommendation of a licensed healthcare professional. High intensity PEMF (pulsed electromagnetic fields) should not be used or viewed as a treatment or cure for any cancer related conditions. The effects of PEMF on this, have not been fully understood and therefore deemed a contraindication.

Metal Implants in the Treatment Area

While many modern orthopedic implants (e.g., titanium screws, plates, or joint replacements) are considered non-magnetic, high-intensity PEMF may still cause localised heating or discomfort. Avoid placing applicators directly over such implants unless you are cleared to do so by a medical professional.

Note: If you are unsure whether any of these contraindications apply to you, always consult with a qualified healthcare provider before using the CELLER8® Max. Your safety is paramount.

Precautions

Keep the device at least 2 meters away from other high-powered electronics (computers, Wi-Fi routers, medical diagnostic equipment) while in use.

High-intensity PEMFs may cause mild dizziness, temporary fatigue, or a metallic taste in sensitive individuals. Discontinue use if symptoms persist.

Consult a healthcare professional before use if you have a serious medical condition or concerns.

Remain well hydrated before and after use. Adequate hydration supports comfort and general wellbeing during sessions.

Not suitable for children under 18 without medical advice and supervision.

Do not use while operating heavy machinery or driving.

Keep away from water and damp environments. Use in a well-ventilated, indoor environment.

Do not operate if the device is damaged.

Use only as directed, and only use accessories approved by the manufacturer.

If you have any concerns or require further information, please contact the CELLER8® team at office@celler8.com. Do not use if you are unsure whether the device is safe for you.

Package Contents

1. CELLER8® Max
2. Touchscreen and control panel
3. Applicator port A
4. Applicator port B
5. Handle
6. Wheels
7. Body mat
8. Large loop
9. Double small loop
10. Applicator carry case
11. Power cord
12. Remote controller



Using the CELLER8[®] Max

Do not power the device on until it is in the location of a session, the applicators are connected and in position and the user is ready for their session to begin. Do not use the CELLER8[®] Max if you have not read this manual in full, especially the contraindications, warnings and safety information. Never touch the applicator ports while the device is on

Opening the protective casing

Always store or relocate the CELLER8[®] Max with the casing closed.

To open the case, apply light pressure to the two clips on the casing. With your fingers on the clips and your thumbs against the casing, pull the clips towards you and push away with your thumbs. Then open the front casing, ensuring the device is standing upright. Once the casing is open and the device is in the location to use, you may connect the applicators.

Connect the applicators

Choose the applicators you wish to apply during the session. You can either connect one applicator or two applicators at the same time.

Put the applicators in position, ready for the session. Then you can locate the applicator ports on the CELLER8[®] Max. Make note of which applicators you connect to applicator port A and applicator port B.

To connect the applicators, push them onto the port connectors and only by the correct colour. It is very important to match the right colour to the right connector, or you could damage the applicators and device. Push them on hard, until you hear them click. If they have not clicked in, you have not pushed them in securely enough.

To remove the applicators, if you want to change them or store them after use, use a finger and thumb to squeeze the buttons on either side of the applicator ports, and then pull the connector to remove it safely. These can then be stored in the applicator box, in a safe location.

Plugging in the power cord

Only once the CELLER8® Max device is situated in the position it will be used in, and the applicators are also safely in place as well as securely connected to the CELLER8® Max, can the power cord be connected.

Please ensure the wall socket is in the 'off' position before it can be connected. Then you can connect the power cord to the CELLER8® Max device, on the front panel, located close to the red power switch. Push the power cord in securely until it will no longer go any further. There will be enough resistance to hold the power cord in by itself when it is secured properly. Then connect the plug end to the wall socket.

Please check everything is in the correct location and when you and the user are ready to start a session can you turn the wall socket to the 'on' position.

Turning the device on

To turn the device on once all the steps above have been followed safely, switch the red switch on the front panel of the CELLER8® Max. The switch should illuminate red and you will see a 'WELCOME' message on the touch screen.

WELCOME



Navigating the menu

Touch the screen to begin. You will see the main  screen, with the CELLER8® Max logo.

Here, from left to right, you can:



Turn the sounder off.



Return home.



Set up a session in settings.



View the protocols.

Pre-Set Protocols


To view the pre-set protocols in preparation for a session, press the  button and you will see the protocols listed in the table.

To use a protocol, simply press on the protocol you wish to use, this will then automatically set up your session and take you to the session page.

These are just a guide, and do not form the basis to every session or guarantee the success of a session. These are listed to guide and take the guesswork out of the initial sessions.



	INTENSITY	FREQUENCY	TIME
RELAX	10%	2	5 MIN
RECOVER	20%	10	10 MIN
ACCELERATE	40%	15	15 MIN
NECK	10%	2	10 MIN
SHOULDER	20%	6	10 MIN
BACK	20%	6	10 MIN
ABDOMEN	30%	4	10 MIN
PELVIC	20%	4	10 MIN
KNEE	30%	8	10 - 15 MIN
THIGH	30%	8	10 - 15 MIN
FOOT	20%	4	10 - 15 MIN
MUSCLE TONE	50%	20	10 - 20 MIN

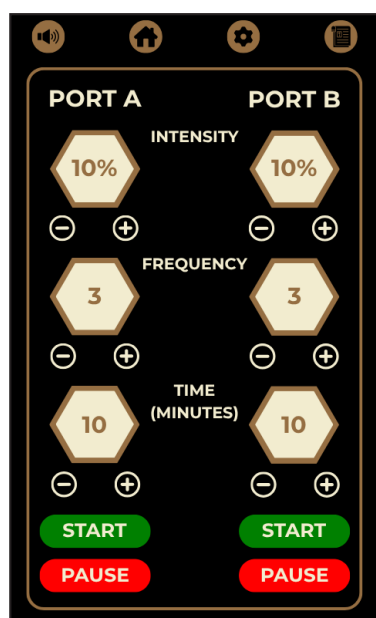

Setting up a session

To set up a session, press the  button.

Here you can change the **intensity**, **frequency** and **time**. Take careful note which applicator you are changing the parameters for. Applicator A at the top, relates to the setting directly below them. These are also clearly marked on the front panel of the CELLER8® Max. The intensity, frequency and time (minutes) can all be increased or decreased, using the  or  buttons.

Starting a session

Before starting any session, the user must remove all jewellery, watches, electronic equipment, credit cards, and metal clips or buttons

Once you have your desired settings and are ready to begin a session. Press  and the CELLER8® Max will begin.

It is **always** encouraged to start at a low intensity setting and slowly increase it throughout the session, for comfort and safety of the user. Never reach an intensity that is uncomfortable for the user during a session. New users should not do more than 5 minutes for their first session. The average session thereafter should be 10 minutes. During the session, the time (minutes) will continue to fall, each minute, until the session has ended.

Pausing a session

To pause a session, press the  button and the device will stop pulsing. You can keep the device paused until you are ready to resume. At which point you can press  again to resume the session.

Ending a session

To end or finish a session early, press the  button and then allow the device to cool for at least 5 minutes before turning the device off with the red power switch. Alternatively, leave the device to end the session once the timer has counted down to zero. The pulsing will stop automatically and then leave the device to cool for at least 5 minutes before switching it off with the red power switch.

Continuous Pulse

Press the  button to activate the continuous pulse mode, which delivers uninterrupted pulses with no pause between frequency cycles.

Continuous Mode should not exceed 20 minutes of use. **Misuse of this mode will void your warranty.**

Turning the sound on/ off

To turn the sound off, press the  button at any time. A line will present through the image of a speaker,  indicating the sound is now off. To turn the sound back on again, just repeat this process until the line has disappeared and the  button is now presented.

Turning the device off

Always leave the device for at least 5 minutes after use to allow the fans to cool the device, before powering down. To turn the device off and to ensure the device is off, press the red switch to the off position and there will be no illumination, to indicate the device is now powered down. When not in use it is best practice to turn the plug socket off and unplug the device from the wall socket completely.

Using the device multiple times per day


Always allow at least 20 minutes between sessions, to avoid overheating. Overuse or using too frequently can result in damage to the device and this is deemed as misuse, and will void your warranty immediately.

Remote Control Function

Please ensure to fit a 3V, coin cell (CR2025) battery to your remote before attempting to use it. This battery is not supplied with the device and must be fitted before first use. You can access the battery compartment at the back of the remote.

Switching On and Linking the Remote

To use your remote, ensure the CELLER8[®] Max device is switched on and at the custom program screen after pressing the  button and by following the instructions previously in this manual. Once at this screen, the remote will already be linked to your controller and with the battery fitted, and will automatically respond to your command.

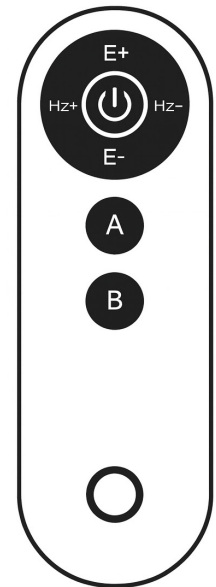

If your remote does not respond while you are on the custom program set-up screen, you will need to manually re-link it to the CELLER8[®] Max. To do this, turn the device on and, as soon as the 'WELCOME' page appears, press and hold the bottom right corner of the touchscreen. While holding this position, press the power button on the remote approximately ten times. Please note that the device will not display any confirmation of connection. To confirm that the process has worked, simply navigate to the custom program screen and test the remote buttons to ensure the device responds.

Assuming your controller is already linked, and the battery is fitted and responds to your commands on the set up screen, you can now use it to set up a session. If the sound is ON , you will always hear a “beep” each time you press a button on the remote. When the sound is OFF , you will not hear a sound but the device will still respond and change the parameters each time.

Selecting Applicator Ports (A & B)

To start, first select an applicator port that you wish to edit or customise by using the (A) or (B) buttons, which correspond with the correct applicator ports on the CELLER8[®] Max. This will ultimately depend on what applicator ports you have already installed applicators to, before turning the device on. If you only have one applicator in port A, select the (A) button and input your desired session before starting a session. If you have both applicator ports with applicators plugged in, then always start with port (A) and then (B) to avoid confusion and ensure you input the correct parameters for your desired session.

Adjusting Frequency (Hz) & Intensity (E – Energy)

Once you have selected the applicator port you want to edit from (A) or (B) on the remote, you can then increase or decrease the Frequency, by using either the (Hz+) or (Hz-) buttons. This will increase or decrease the Frequencies in Hz by 1 increments each time. The (Hz+) button will increase by 1Hz and the (Hz-) button will decrease by 1Hz each time it is pressed.

Alternatively, to increase or decrease the intensity each time, you will need to use the (E+) or (E-) buttons on the remote. 'E' stands for, and represents 'Energy' from the device. Therefore when pressing the (E+) button, it will increase by 10% intensity or decrease by 10% intensity when pressing the (E-) button.

Starting or Pausing a Session

Once you have reached the desired setting for your session from both or just one applicator port. You are now ready to begin by pressing the  button. PLEASE NOTE: You cannot start/ pause both ports at the same time with this button. You must start or pause each port individually by selecting which port from (A) or (B) from the controller and then pressing  each time. You can start or pause an applicator port at any time during your session individually.

You can also use the remote during a session as this makes it easier for the user increase or decrease the intensity or frequency, mid-session, without the need to manually make these changes from the main CELLER8[®] Max controller. This allows for easier application and less need for movement from the user during the session when they wish to make a change to the parameters. PLEASE NOTE: always start at the lower intensity settings and work your way up and not the other way around, as the “low and slow” approach to high intensity PEMF sessions.

Time Settings and Cooling Periods

You are not able to change the time parameters using the remote and it will only run for the default time set on the device, unless you manually change this on the screen of the CELLER8[®] Max before starting your session. Once you have finished your session and the time on the device has run to 0, the session will end. Please allow at least 5 minutes of cooling time before manually turning the device OFF on the red power switch at the front of the device and then at the wall outlet. If you wish to use the device again for another user, please leave at least 20 minutes of cooling time before starting another session.

Applying The Applicators




All applicators must be at least 2 meters away from any other electronic equipment. The users must remove all jewellery, watches, electronic equipment, credit cards, and metal clips or buttons, at all times, during the session.

The **CELLER8® Max body mat** should be laid flat and used on a hard surface. At least 2 meters from any other electronic equipment or monitoring equipment. There is a cover provided for comfort and cleanliness. The covers and mat can be cleaned by following the cleaning guide in this manual.

The large single loop applicator can be applied to a larger area or used with the body mat at the same time. To intensify a session, the user can hold the loop whilst laying on the body mat. This will create a sandwich effect and maximise the power of the session for the user in-between both applicators.

The applicators are not marked North or South, but usually the user will feel the pulses stronger when the polarity of the applicators are synchronised. You can use the large loop applicator with the padel cover on or with just the loop covers on, depending on preference and application.

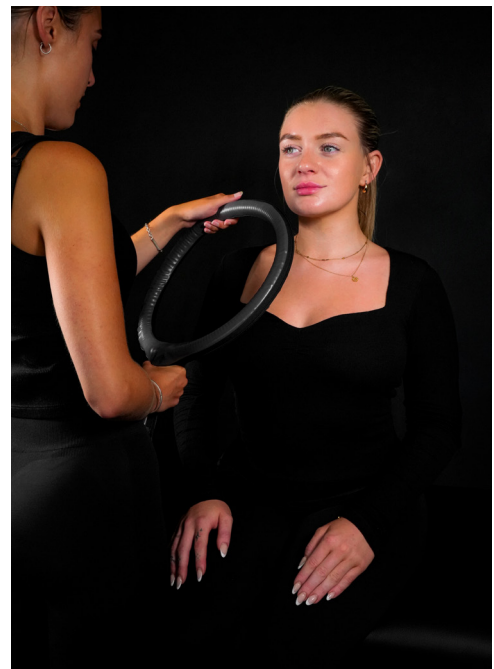

Continued on the next page.

The smaller double loop applicators can be used on their own or at the same time as the body mat. The smaller double loop applicators can be secured together using the small clips at the top of both loops to create a sandwich effect either side of a limb or joint.

Alternatively the loops can be separated and used in two different areas simultaneously. Again, the applicators are not marked north or south, but the user can usually feel one orientation will be stronger than the other, indicating the polarities are matched.



CELLER8[®] Max Technical Specifications

Specification	Details
Device Type	High-Intensity PEMF (Pulsed Electromagnetic Field) System
PEMF Intensity Range	Mats Intensity Range: Approx. 800 - 8,000 Gauss Rings Intensity Range: Approx. 2,400 - 24,000 Gauss
Frequency Range	1 – 20 Hz
Pulse Type	Rapid rise and fall, high slew rate pulses
Application Modes	15 Suggested Protocols, Program is Customisable
Session Duration	1 – 20 minutes per cycle
Applicators Included	1x Full-body mat, 1x Large Loop, 1x Dual Small Loop
Full-Body Mat Dimensions	Approx. 132 × 52 cm
Device Dimensions	Approx. 58 × 24 × 34 cm
Weight	Approx. 20 kg
Power Supply	Standard AC Power (region-specific) 220-240V 110V or 220-240V available (region specific)
User Interface	Digital controller with touch screen
Construction	Solid-state cabinet housing
Use Case	Personal or professional wellness use
Classification	Non-medical device, CE
Power Rating	566W - 850W

Maintenance & Cleaning

Proper maintenance and regular cleaning of your CELLER8® Max device and its accessories will help ensure optimal performance, user safety, and product longevity. Please follow the guidelines below carefully.

General Maintenance

Inspect regularly: Before each use, inspect the device, applicators, and cables for signs of damage, wear, or loose connections. Do not use the device if any components appear compromised or damaged.

Avoid overuse: Allow the device to cool for at least 20 minutes between sessions to avoid overheating and ensure internal components are not stressed.

Cable care: Avoid sharp bends or pulling on cables. Always disconnect cables by gripping the plug, not the cord.

Storage: When not in use, store the device and all accessories in a dry, dust-free area away from direct sunlight, moisture, or extreme temperatures.

Transportation: If transporting the device, use the original packaging or a padded case to protect the unit from impact and vibration.

Cleaning Instructions

Unplug first: Always disconnect the device from the power source before cleaning. Never attempt to clean the unit while it is in operation.

External surfaces: Wipe the outer surfaces of the control unit and applicators with a soft, damp cloth. For more thorough cleaning, use a cloth lightly moistened with a mild soap solution.

Disinfection: For shared or professional settings, you may use non-abrasive, alcohol-free disinfectant wipes that are safe for electronic devices. Avoid excessive moisture near buttons, display, or cable ports.

Do not submerge: Never immerse the device or any of its components in water or cleaning liquids.

Avoid harsh chemicals: Do not use bleach, ammonia-based cleaners, or abrasive pads, as these can damage the surfaces and reduce product lifespan.

Servicing

Although the CELLER8® Max contains no user-serviceable internal parts, annual servicing is encouraged by a CELLER8® authorised representative. This ensures there are no warning signs of wear and tear that could result in any potential hazards. If the device is not functioning correctly, please contact CELLER8® support for troubleshooting or to arrange an authorised service.

Unauthorised repairs or tampering will void the warranty immediately.

Important: Regular care ensures your CELLER8® Max continues to operate effectively and safely. Always follow these maintenance guidelines to maintain your warranty and product performance.

Troubleshooting

The CELLER8® Max is built to deliver reliable performance. However, if you experience any issues during use, refer to the following troubleshooting guide before contacting support.

Issue	Possible Cause	Suggested Solution
Device does not power on	- Power cable not connected properly - Power source is not working - Internal fuse may be blown	- Ensure the power cord is fully inserted into the device and outlet - Try a different wall socket - Contact CELLER8® support if issue persists
No PEMF output	- Applicator not connected correctly - Program not started - Low intensity setting	- Check all cable connections - Ensure you have selected and started a program - Increase intensity setting if appropriate

Issue	Possible Cause	Suggested Solution
Intermittent or inconsistent signal	- Loose cable connection - Cable or applicator fault	- Disconnect and reconnect all applicator cables - Try a different applicator if available - Contact support for replacement parts if needed
Unusual noise or vibration	- Normal operation sounds - Internal fan or pulse switch	- Some clicking or pulsing sounds are normal - If noise is excessive or accompanied by performance issues, stop use and contact support
Device overheating	- Extended use without breaks - Blocked ventilation - High ambient temperature	- Allow the device to cool between sessions - Ensure proper airflow around the unit - Avoid using in hot or enclosed spaces
Display or buttons unresponsive	- Software glitch - Moisture or debris on interface	- Power off and restart the device - Wipe display area with a soft dry cloth - If unresponsive, contact support

Still having trouble?

If you have tried the suggestions above and the issue is unresolved, please contact CELLER8® Customer Support with the following information:

- Device serial number (found on the back or underside)
- Description of the issue
- Any recent changes in usage or environment

CELLER8® Max customer support:

Email: office@celler8.com

Website: www.celler8.com

Telephone: +44 (0) 1795 379044

Monday – Friday, 9:00 – 16:30 (standard UK time zone GMT)

CELLER8® HQ: NewMed Limited

Hop House, Woodgate Lane, Borden, Sittingbourne, Kent, ME9 8JX,
United Kingdom.

Warranty & Support

We stand behind the quality and reliability of our products. The CELLER8® Max comes with a standard limited warranty to ensure your satisfaction and peace of mind.

Warranty Coverage

CELLER8® warrants the CELLER8® Max device and its original accessories against defects in materials and workmanship under normal use for a period of 2 years from the date of purchase.

This warranty covers:

- Manufacturing defects in the main device unit
- Faulty controller or display functions
- Defective cables or connectors (under normal usage)
- Non-functioning applicators (mat or loops)

This warranty does not cover:

- Damage caused by misuse, accidents, or unauthorized repairs
- Normal wear and tear of surfaces, cables, or connectors
- Cosmetic damage (e.g., scratches, dents)
- Damage due to improper storage, transportation, or handling

Warranty Claim Process

To make a warranty claim, please contact CELLER8® Customer Support and provide the following information:

- Proof of purchase (e.g., order confirmation or receipt)
- Device serial number (located on the unit)
- Description of the issue and any troubleshooting steps taken
- Photos or videos if applicable

If the issue falls under warranty, we will provide instructions for return, replacement, or repair as applicable.

Customer Support

If you have any questions, technical issues, or need assistance with your CELLER8® Max, please reach out to our support team:

Email: office@celler8.com

Website: www.celler8.com

Telephone: +44 (0) 1795 379044

Monday – Friday, 9:00 – 16:30 (UK Time)

CELLER8® is a brand name, owned and distributed worldwide by NewMed Limited.

CELLER8® HQ:

Hop House,

Woodgate Lane,

Borden, Sittingbourne,

Kent, ME9 8JX,

United Kingdom.

